


myVIBS - OPERATING INSTRUCTIONS

CONTENT

- I. Components of *myVIBS*
- II. Safety Instructions
- III. Starting/Operation of *myVIBS*
- IV. Usage
- V. Cleaning and Maintenance
- VI. Inserting/Replacing the Battery
- VII. Warranty
- VIII. Protection of the environment

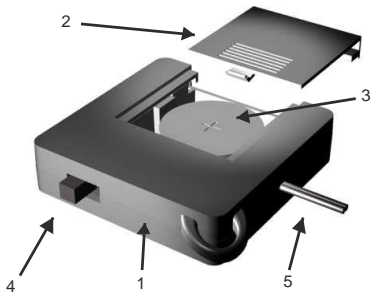

I. COMPONENTS OF *myVIBS*

- 1 Casing
- 2 Battery Cover
- 3 Battery
- 4 Switch off/on
- 5 Cable
- 6 Mouthpiece Vibrator
- 7 Transport bag (not shown)
- 8 Lanyard (not shown)

II. SAFETY INSTRUCTIONS

- Don't expose the unit to temperatures below -5°C/23°F or above +40°C/104°F for longer periods of time.
- Check that the battery is in right position (+ is in upper position).


- Don't squeeze the cable (5).
- Don't try to repair the unit yourself.
- Never submerge any component of the unit into a liquid as this will cause a short circuit and destroy the electronics.

III. STARTING /OPERATION OF *myVIBS*

- Before starting please insert the supplied coin battery 3 according to VI into the battery slot of the casing 1.
- Battery life should be about 6 hours. The number of off/on cycles however may change the lifetime of the battery (3).
- It is recommended to use *myVIBS* for about 8-10 minutes/day.
- *myVIBS* never switches off automatically: it is necessary to switch it off/on manually with the switch 4 on the casing 1.

IV. USAGE

- Place the mouthpiece (6) between your teeth (c, d) and lips (a, b).
- Try to place the mouthpiece in a comfortable position where it does not strain the muscles around the lips and does not touch the teeth, i.e. try to maintain a separation zone (e), eventually by slightly extruding the lips.
- Switch the button on/off.

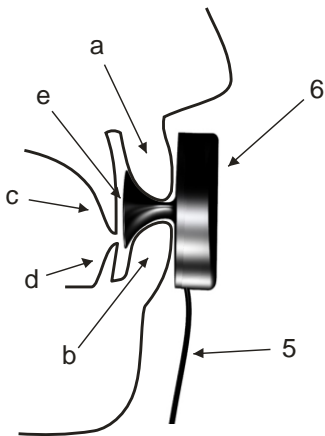

- a Upper lip
- b Lower lip
- c Upper front teeth
- d Lower front teeth
- e Separation space

- ***For relaxation:***

Keep your lips closed and maintain some separation space between the teeth and the lips (eventually extrude the lips to a position as if you want to kiss).

- ***For training:***

Try to do exercises with your lips (e.g. smile / relax cycles) while *myVIBS* is operating.

- It is recommended to use *myVIBS* daily; decide yourself what you need: relaxation or training.

V. CLEANING AND MAINTENANCE

- You can clean the mouthpiece with appropriate products, e.g. using towels with liquid soap or refreshing towels or disinfection towels.
- The casing can be cleaned with anti-static cloth.
- Never submerge the mouthpiece into water or any other liquid, because this will destroy the vibrator electronics.
- To protect *myVIBS* from mechanical damages use the bag provided in the supplied presentation box.
- Use the supplied presentation box for storage when not using the unit for longer periods of time.

VI. INSERTING/REPLACING THE BATTERY

- To replace the battery 3 open the battery cover 2 (press your finger on the marked zone of the battery cover 2 and push the cover carefully along the surface to the outside of the casing).
- Remove the old battery with your fingernail/a ballpen/a screwdriver.
- Insert the new battery (battery type: CR2450) making sure the + pole is facing up and press on the battery so that it is fixed inside the box between the electrical contacts.
- Finally close de battery cover 2 carefully and let the snapper snap in while slightly pressing on the marked zone.

VII. WARRANTY

- *myVIBS* has 1 year of warranty on the material that has been provably defect from the factory, starting the day of sales.
- For requesting the warranty it is necessary to send a message with the details of the defects to: costumer@myvibs.com
- Verbal o telephone communications are not sufficient.
- The warranty is carried out exclusively for the original buyer.
- Battery is excluded from the warranty.
- The warranty is invalid in the follow cases:
 - The damage has been caused by inappropriate handing.
 - Repairs have been made by an unauthorized workshop.

VIII. PROTECTION OF THE ENVIROMMENT

- Used batteries must be disposed of through collection bins marked for this purpose.
- Never put the battery in the trash.
- Dispose of the *myVIBS* apparatus according to European Directive 2002/96/EC on Waste Electrical Devices.
- For more information contact the municipal authority responsible for waste.





myVIBS
Never forget me!

INEXPIM S.L.
Antonia Merce, 2
28009 Madrid/Spain

www.myvibs.com